

ON THE MEASUREMENT OF EXTENSIONAL VISCOSITY OF POLYMER MELTS USING SER UNIVERSAL TESTING PLATFORM

Filip P.* and Svrčinová P.
*Author for correspondence
Institute of Hydrodynamics,
Acad.Sci.Czech Rep.,
Prague 6, 166 12,
Czech Republic,
E-mail: filip@ih.cas.cz

ABSTRACT

Measurement of elongational viscosity of polymer materials seems to be much more complicated in comparison with measurement of shear viscosity. The problems concern preparation of the polymer samples (i.e. their homogeneity and geometrical accuracy), fixing of both ends, influence of gravitational force, etc. Usage of the individual devices for measurement of elongational viscosity differs, among other things, by an aspect ratio of the prepared samples. In contrast to 'classical' devices developed in 80s' using oblong samples, a recently proposed SER Universal Testing Platform enables to apply rectangular samples of moderate aspect ratio. However, in context with this ratio, there has been recently discussed [1] a possible appearance of planar elongation during measurement of uniaxial elongational viscosity. The aim of this contribution is to document - for material used (LDPE Escorene, Exxon, USA) - the invariance of the elongational viscosity measured with respect to sample widths and thicknesses (within sufficiently broad regions). Temperature was set to 180°C.

INTRODUCTION

Viscosity represents one of the main characteristics of polymer melts, therefore rheological measurements focus on determination of this material property in different flow geometries. In polymer processing polymeric materials usually subject to two types of deformation: shear flow and extensional flow. This results in two corresponding types of viscosities: shear and extensional.

Shear flow is the most common type of flow discussed in rheology. Shear flow may be produced in many ways and determination of shear viscosity is relatively well documented. In some polymer processing operations as e.g. fibre spinning or injection moulding, there occurs a combined flow exerted by both shear and elongational forces. In this case mere shear

experiments cannot describe real polymer flow and a proper analysis of extensional flow is inevitable.

Extensional flow measurements are very sensitive to molecular structure, crystallinity and the macrostructure effect (mainly influenced by long-chain branching) of the polymeric system being tested. The measurements face certain difficulties as necking and rupture or sagging during heating a sample. Often it is not possible to reach a steady stress before the sample ruptures or irregularly deforms. A couple of methods have been tried to solve these problems and generate purely extensional flow. There exist several instruments such as a rotating clamp or extensional rheometer with a translating clamp and vertical (or horizontal) buoyancy control bath, proposed by the leading researchers in this branch, e.g. Meissner (1969), Münstedt (1975,1979), Meissner et al. (1982).

In recent years there has been developed a new instrument, known as a SER (Sentmanat Extensional Rheometer) Universal Testing Platform. This detachable unit is used with commercially available rotational rheometer systems (TA Instruments, Physica series MCR, etc.). The SER unit consists of paired master and slave wind-up drums on bearings housed and mechanically coupled via intermeshing gears. Rotation of a drive shaft results in a rotation of the affixed master drum and rotation of the slave one with equal angular velocity but in the opposite direction. This causes the ends of the sample to wind up onto the drums, and thus stretches a sample over an unsupported length. The ends of the sample are attached to the drums by two securing clamps as documented in Fig.1, and described in detail by Sentmanat (2003,2004), Sentmanat & Hatzikiriakos (2004).

The present contribution aims at the measurements of transient uniaxial elongational viscosity of polyethylene material using this device including a discussion on advantages and disadvantages of this experimental method.

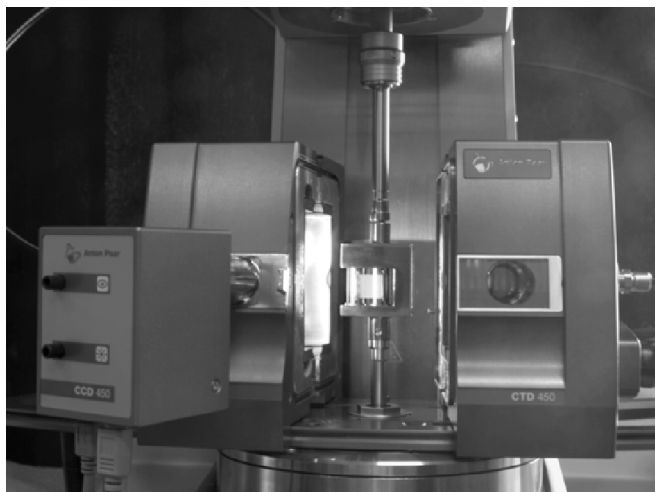


Figure 1 Experimental set-up (SER Universal Testing Platform surrounded by both halves of Convection Temperature Control Device CTD450) with a sample at room temperature and a detail of SER Universal Testing Platform with an individual clamp

EXPERIMENTAL

Experimental device

Measurements were carried out using the SER Universal Testing Platform from the Xpansion Instruments; the model SER-HV-P01 was applied with a Physica MCR501 rotational rheometer host system. For measurement of rheological properties of polymer melts there was used a convection heated measuring chamber CTD450 equipped with a camera system. The images and videos were recorded and displayed directly in the Rheoplus software supporting Physica MCR501.

Material

The low density polyethylene (LDPE) Escorene used in the experiments is produced by Exxon Mobil (basic characteristics in Tab.1). Sample specimens of LDPE were prepared by compression moulding in a preheated press. The final dimensions of the solid polymer rectangular specimens were approximately 12.7 mm (1/2") in width and 18 mm in length, with thickness in the range of 0.6-0.9 mm.

Table 1 The basic characteristics of the used material.

Material	MFI [g/10min]	Density [g/cm ³]	M_w [g/mol]	M_n [g/mol]
LDPE Escorene (Exxon, USA)	0.33	0.922	366,300	30,280

Measurements

Transient uniaxial extensional viscosity was measured at temperature 180°C. For characterising polymer melt there were conducted start-up experiments with constant strain rates (within the range 0.0316-10 s⁻¹); the results are shown in Fig.2. If a sample rupture appears if one revolution of the drums is not yet attained, then there is no problem with determination of the magnitudes measured. A little complication occurs if a sample has to surmount the individual clamps (more than one revolution). This phenomenon is accompanied by a diversion (steep peak) of an experimental curve (very roughly resembling δ function) from the expected course and as such can be transformed to true values by continuous connection of the preceding and subsequent parts of the curve investigated. However there is no compensation if this ' δ -clamp behaviour' is exhibited at a maximum as demonstrated in Fig.2 for extension strain rate 0.0316 s⁻¹. To eliminate this case there is a possibility to change a sample thickness.

Fig.3 depicts behaviour of extensional viscosity vs. Hencky strain for six different strain rates within the range 0.0316-10s⁻¹.

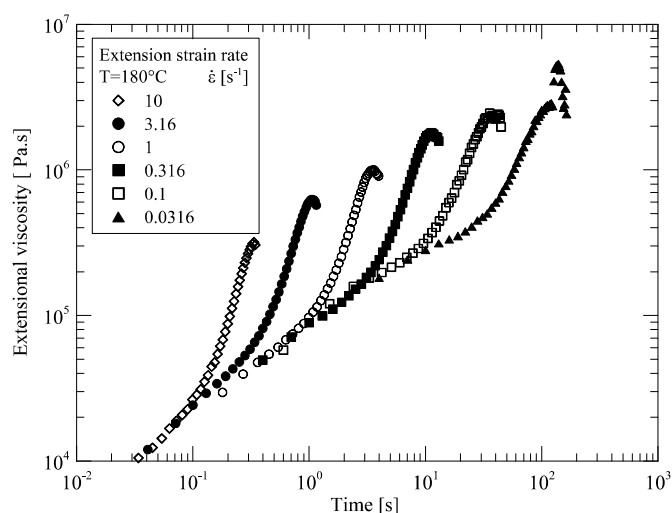


Figure 2 Behaviour of transient uniaxial extensional viscosity

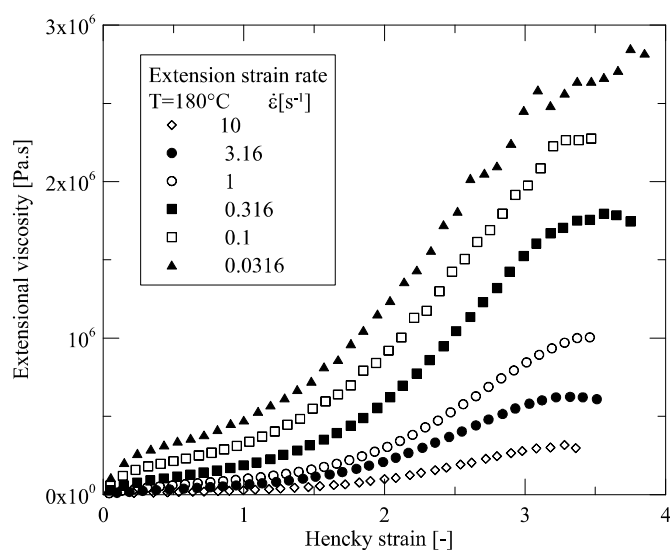


Figure 3 Transient extensional viscosity as a function of Hencky strain at six different strain rates.

DISCUSSION

For proper measurement of elongational viscosity using the SER unit there should be fulfilled the following assumptions (or at least almost fulfilled):

- 1) even distribution of temperature throughout a whole sample during measurement;
- 2) a sample should be stretched along its full width;
- 3) prevention from degradability.

ad 1) Due to its geometrical arrangement (including the air inlets and outlets) the oven CTD450 exhibits a uniformly distributed temperature field even in the presence of the asymmetric SER unit. This results in a uniform heating of a whole sample and a sample thickness only influences a time period for which the sample should be exposed to reach the pre-set temperature.

ad 2) A rectangular shape (see Fig.4) is an optimal one for starting elongational measurement as this ensures proper stretching of the whole sample along its width just from the beginning. To achieve this shape is not so trivial as gravitational force acting on an unsupported part of a sample during heating causes its sagging. For observation of this phenomenon the oven CTD450 is equipped with a camera enabling direct checking of the whole process including its recording. This provides a priori an elimination of the measurements for which a location of the sample measured is unfavourable and should result in an incorrect data. Fig.5 illustrates an application of the camera system for the selected points (however the recording is continuous during the whole measurement).

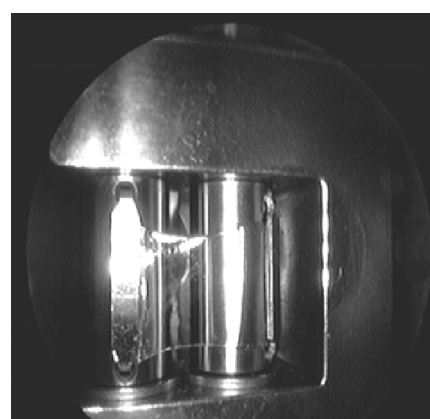
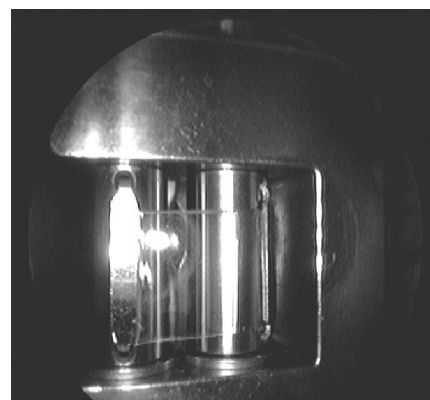


Figure 4 Shape of the sample:
top - optimal shape, bottom - sagging of the sample

ad 3) Time required for proper heating of the whole sample is not limited only by a possible sagging as discussed above but by material degradation as well. This can be partially checked using a camera (colour changes). There are two possibilities how to accelerate the measurements and hence to suppress material degradation:

a) Usage of a thinner sample should shorten a time period for uniform heating to a required temperature. However due to other factors it is recommended not to reduce a sample thickness considerably.

b) Possible elimination of the clamps used for fixing the ends of a sample measured exhibits two favourable effects - both this accelerates measurement itself (and thus postpones possible onset of degradation) and substantially eliminates the 'δ-clamp behaviour' discussed above. This removal of the clamps is not applicable for all materials, nevertheless for LDPE Escorene that was measured the sticking of the ends of samples directly to pre-heated drums is admissible as proved in Fig.6.

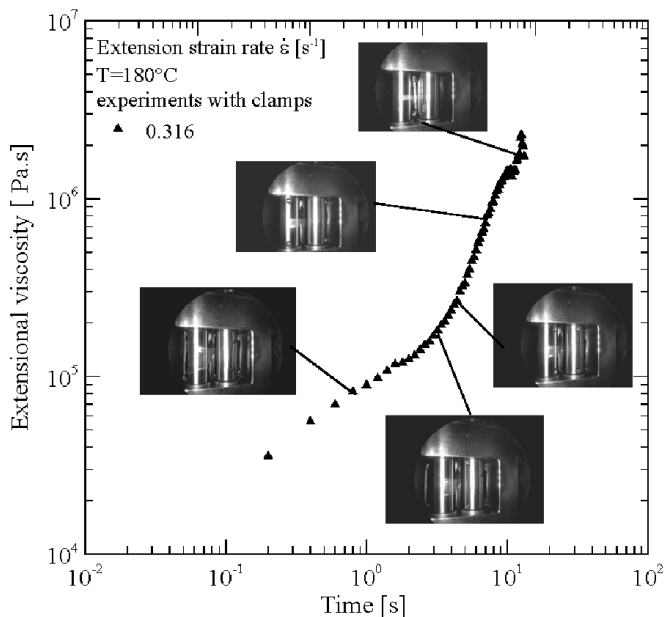


Figure 5 Geometric shapes of the sample for a selected value of extension strain rate

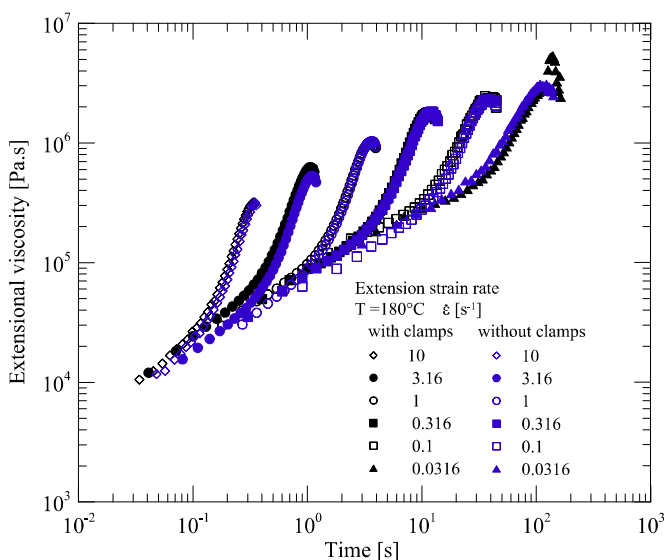


Figure 6 Comparison of uniaxial extensional viscosity measured with and without the clamps

CONCLUSION

The experimental set-up discussed above provides more responsibility for proper evaluation of extensional measurements due to visual checking, including a movie record allocating corresponding pictures to the individual measured points and eliminating unexpected material degradation.

There was shown that the experiments can be substantially accelerated if the clamps are not used.

On the other hand measurement of the transient elongational viscosity using the SER Universal Testing Platform has some limitations and other measurement instruments [1-4] can provide wider possibilities in some respects as e.g. better temperature control.

However for majority of standard materials the results measured using various instruments are fully comparable.

ACKNOWLEDGEMENT

The authors wish to acknowledge GA CR for the financial support of Grant No.103/08/1307 and the Institutional Research Plan No. AV0Z20600510.

REFERENCES

- [1] Yu, K., Marin, J.M.R., Rasmussen, H.K. and Hassager, O., Modeling of Sentmanat Extensional Rheometer, *5th Annual European Rheology Conference*, 9-11 April, 2009, Cardiff, U.K., CD ROM.
- [2] Meissner, J., Rheometer zur Untersuchung der deformations-mechanischen Eigenschaften von Kunststoff-Schmelzen unter definiert Zugbeanspruchung, *Rheol. Acta*, Vol.8, 1969, pp.78-88.
- [3] Meissner, J., Stephenson, S.E., Demarmels, A. and Portman, P., Multiaxial elongational flows of polymer melts - classification and experimental realization, *J. Non-Newtonian Fluid Mech.*, Vol.11, 1982, pp.221-237.
- [4] Münstedt, H., Viscoelasticity of polystyrene melts in tensile creep experiments, *Rheol. Acta*, Vol.14, 1975, pp.1077-1088.
- [5] Münstedt, H., New universal extensional rheometer for polymer melts - measurements on a polystyrene sample, *J. Rheol.*, Vol.23, 1979, pp.421-436.
- [6] Sentmanat, M.L., Dual wind-up extensional rheometer, *US Patent* No 6,578,413 (2003).
- [7] Sentmanat, M.L., Miniature universal testing platform: from extensional melt rheology to solid-state deformation behaviour, *Rheol. Acta*, Vol.43, 2004, pp.657-669.
- [8] Sentmanat, M.L. and Hatzikiriakos, S.G., Mechanism of gross melt fracture elimination in the extrusion of polyethylenes in the presence of boron nitride, *Rheol. Acta*, Vol.43, 2004, pp.624-633.